




## Set Your Own Rhythm for Recovery

InTandem uses a combination of music, foot sensors, and software to help improve walking.

### Get Started

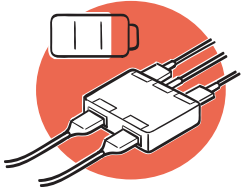

Charge the system for at least 1 hour.



Connect to WiFi or use the preloaded cellular plan to get the latest music and updates.



Get familiar with the device and user manual.

### Your First Walk



Choose a safe place to walk that's flat and level.



Turn on the touchscreen device and attach the foot sensors.




When the music starts, all you need to do is walk!

### Questions? Contact Us

**(844) 400-9255**

[intandemrx.com/support](https://intandemrx.com/support)

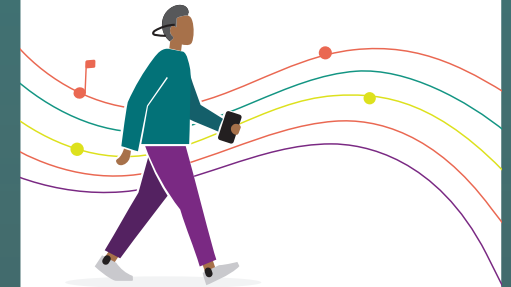

Manufactured by:  
MedRhythms™, Inc.  
183 Middle Street, Suite 300  
Portland, ME 04101

Review the included Instructions for Use for detailed instructions, warnings, and precautions. IFU002, Rev. B, DEC2023.

© 2023 MedRhythms, Inc. All rights reserved. MedRhythms and InTandem are trademarks of MedRhythms, Inc.